



Your partner
in your everyday
surgery

www.neuro4vets.be



NEURO & SPINAL SYSTEMS

OVERVIEW



LeiMESH

3D Dynamic Mesh System, Titanium

A versatile and customizable solution for cranial and maxillofacial bone reconstruction, LeiMESH can be easily cut, shaped, and contoured to match individual anatomy, providing excellent structural stability.



C-LOX

Cervical Intervertebral Fusion & Distraction Implant, Titanium

An intervertebral fusion implant for the cervical spine in dogs, designed to treat Wobbler Syndrome by replacing the spinal disc and enabling rapid bone ingrowth through its open-porous titanium structure.



LeiCAGE

Intervertebral Cages, Titanium

An intervertebral titanium cage system designed to create distraction, support spinal fusion, and ensure reliable osseointegration with excellent biocompatibility and imaging clarity.






LeiPED

Pedicle Screw System, Titanium

A titanium pedicle screw system for spinal stabilization in dogs, designed to treat DLSS, vertebral instability, and spinal trauma through bone-near rod placement and an asymmetric, side-loading tulip head that improves biomechanical stability.

LeiMESH

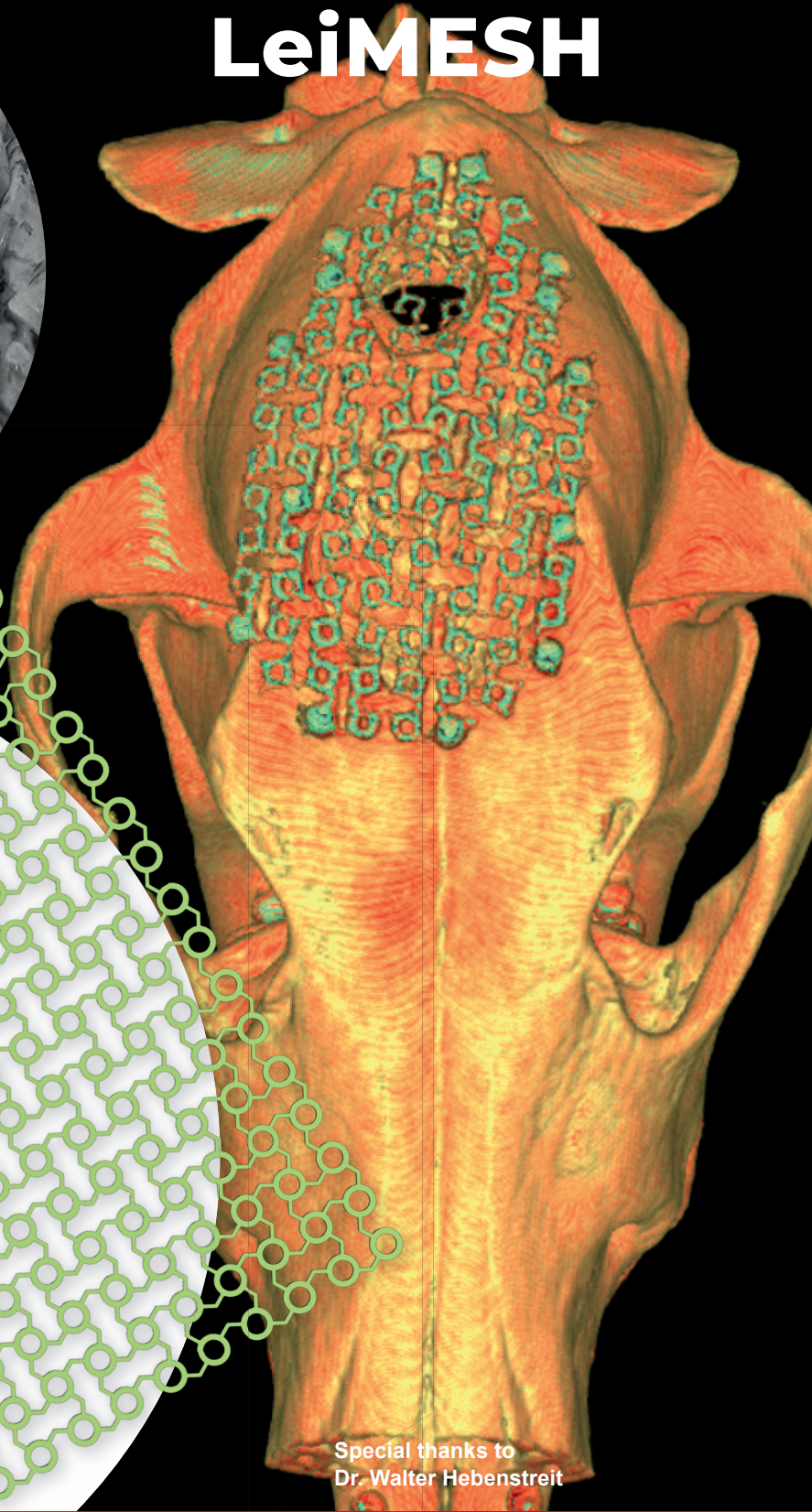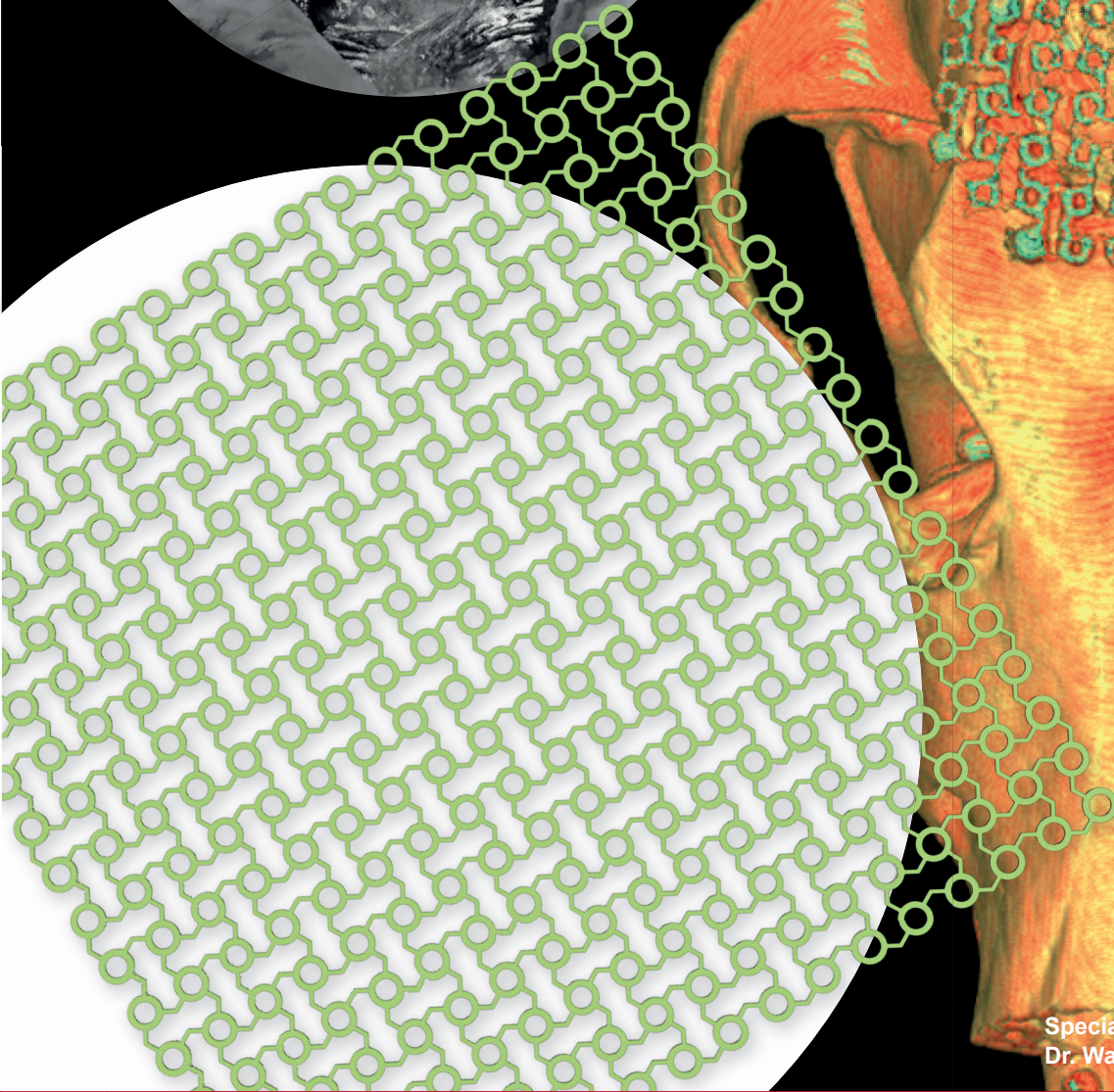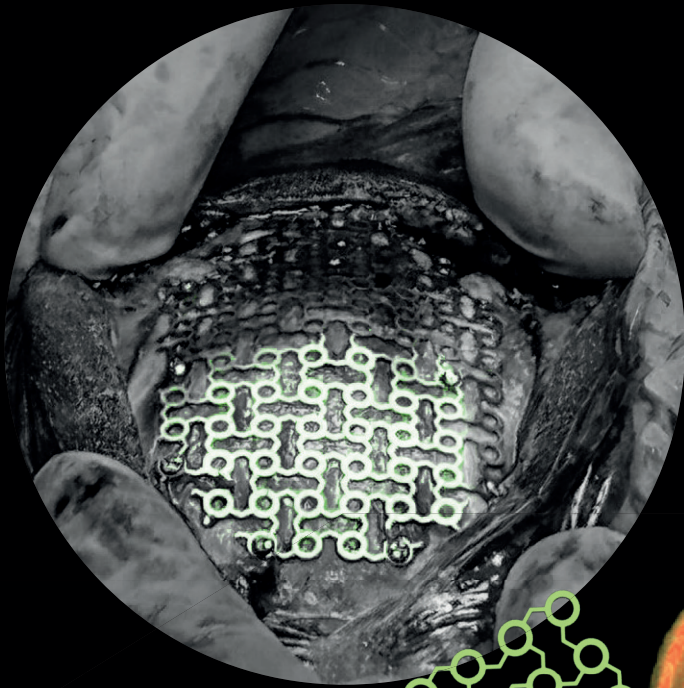

Special thanks to
Dr. Walter Hebenstreit

DYNAMIC MESH SYSTEM

R RITA
LEIBINGER
MEDICAL


LeiMESH

Cuttable Mesh Plate System




Three-Dimensional Adaptability

The Titanium LeiMESH features a unique geometric pattern design. This creates a flexible structure that can be easily bent and shaped, which is essential for creating a custom fit during surgery. It becomes more rigid when formed into a three dimensional configuration, providing stable support.

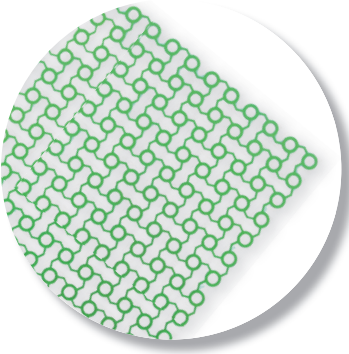

Titanium for Best Biocompatibility

Made of medical grade Titanium, LeiMESH is highly biocompatible, strong yet lightweight, corrosion-resistant, and promotes bone integration. Its flexibility makes it suitable for custom shaping and its non-magnetic properties makes it safe for MRI scans.

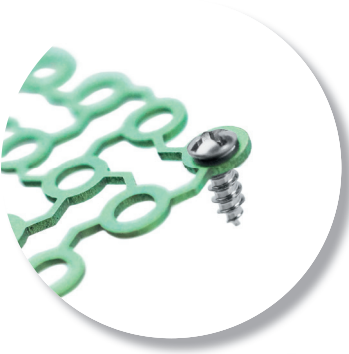

Perfect Balance of Strength and Flexibility

The thickness of LeiMESH provides sufficient structural support while allowing flexibility to customize the mesh to fit complex contours.



Various Applications

LeiMESH is highly versatile and can be used in various medical applications such as cranial fracture repair, bone fracture management, bone cement augmentation, and maxillofacial reconstruction surgeries.



Low Profile Cross Drive Titanium Screws

Made of Titanium to ensure seamless biocompatibility, LeiMESH Screws feature a low-profile head for minimal tissue irritation. Compatible with the LeiMESH Screwdriver, the LeiMESH Cross Drive screws provide a secure fixation, fitting perfectly into the screw holes of the Titanium Mesh to meet diverse anatomical and surgical requirements. Additionally, they offer a self-holding capability similar to LeiStar screws.

LeiMESH 3D Dynamic Mesh Set

LeiMESH 3D Dynamic Mesh Set


Contains:

- 1 LeiMESH 3D Dynamic Titanium Mesh
- 5 of each LeiMESH Titanium Screws
(Ø 1.6 mm, Lengths: 4 - 8 mm // 25 total)
- 1 LeiMESH Screwdriver Handle
- 1 LeiMESH Screwdriver Shaft
- 1 of each Micro Twist Bone Drill
(Ø 1.0 mm, Drill Lengths: 3 mm, 5 mm, 8 mm)
- Universal Wire Scissors
- 1 Tray for LeiMESH Set with Lid

141-9090-00

Tray without contents

141-9090-10

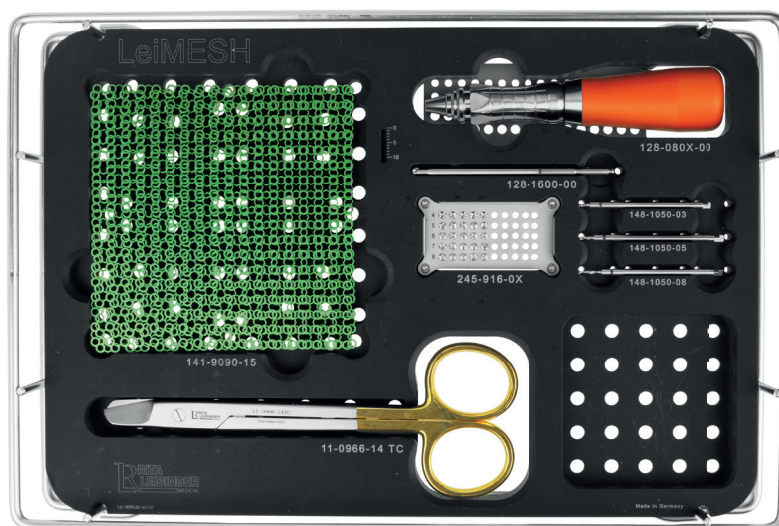

Bending Pliers with Tray



Contains:

- 1 LeiMESH Bending Pliers
- 1 Tray (without lid)

141-9091-00

Tray without lid and contents

141-9091-10

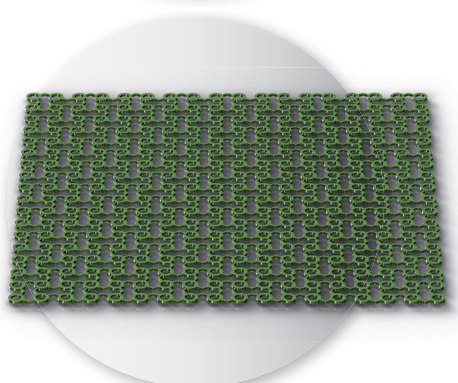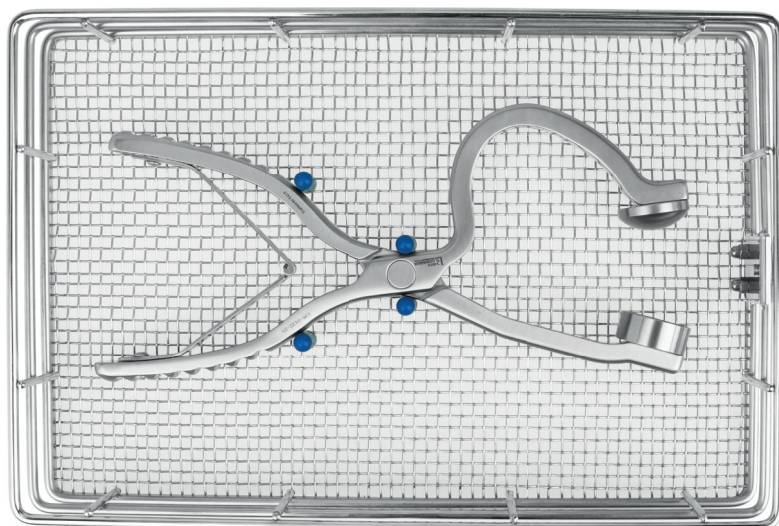

LeiMESH Cuttable Titanium Mesh System

Composed of high-quality titanium, LeiMESH offers surgeons the flexibility to shape implants during surgery.

LeiMESH Implants provide a unique solution for treating cranial or maxillofacial fractures after trauma, reconstruction after tumors, as well as mandibular cases.

LeiMESH Implants are ideal for osteosynthesis in non-load-bearing bones.

LeiMESH Implants are used to augment bone cement in joint replacement procedures.

LeiMESH Implants offer the dual benefits of exceptional biocompatibility and superior corrosion resistance.


LeiMESH Implants are suitable for permanent implantation in animals due to their long-term stability.

LeiMESH Cuttable Mesh Titanium

3D Dynamic Cuttable Mesh Titanium


0.6 mm thick, 90 mm x 90 mm

141-9090-15

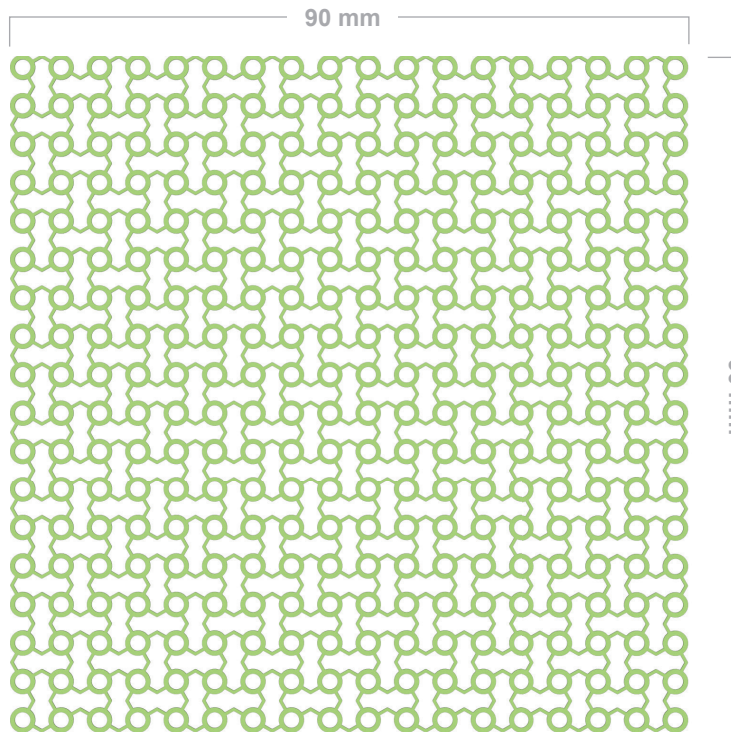

LeiMESH Screw Titanium

Ø 1.6 mm, low profile, Cross Drive

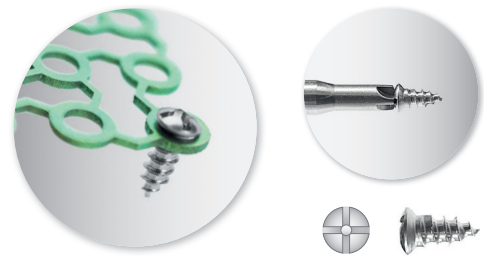

Product Code	Length (mm)
245-916-04	4
245-916-05	5
245-916-06	6
245-916-07	7
245-916-08	8

LeiMESH Screwdriver Handle

with dental connector/coupler

128-0803-00



LeiMESH Micro Twist Bone Drill

Ø 1.0 mm, with dental connector/coupler



Product Code	Drill Length
148-1050-03	3 mm
148-1050-05	5 mm
148-1050-08	8 mm

Universal Wire Scissors

Stainless Steel

11-0966-14TC



LeiMESH Screwdriver Shaft

Cross Drive, with dental connector/coupler

128-1600-00



LeiMESH Bending Pliers

Stainless Steel

128-0122-25



C-LOX[®]
Makes the lame walk again.




Special thanks to
Prof. Björn Meij

C-LOX SPINAL, locking


RITA
LEIBINGER
MEDICAL

C-LOX Cervical Intervertebral Fusion Implant

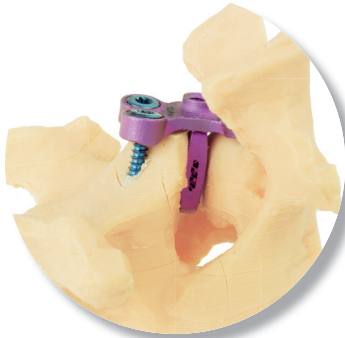

Versatility and Stability

The cages come in different sizes to fit different breeds and spinal segments. Developed to combine the benefits of an intervertebral spacer, rigid titanium plate, and locking screw fixation, the C-LOX Cages not only preserve the natural anatomic profile of the spine, but also provide great stability and convenience in placement and fixation.

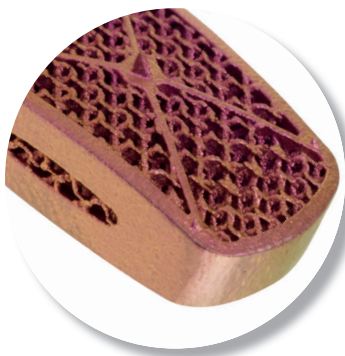

Improved Osseointegration

The dodecahedron construct of the cage provides a larger surface area for bone growth, which improves osseointegration and cell attachment, therefore promoting bone growth and reducing implant failure rates.



Exceptional Biocompatibility & Strength

With minimal adverse reactions or rejections, Titanium reduces the risk of infection and inflammation. Its hypoallergenic properties make it a safe choice, reducing the likelihood of allergic responses. The unique Titanium Alloy also makes our implants and screws very durable but lightweight.



LeiStar Screws

The C-LOX Locking screws have LeiStar screwheads, which allows for best tightening torque. Additionally, the screws provide a self-holding feature when used with a RITA LEIBINGER screwdriver shaft, keeping the screws securely on the driver.



C-LOX Surgical Protocol

C-LOX

Surgical Protocol

Surgical application of a new anchored intervertebral spacer (C-LOX) for the treatment of Canine Cervical Spondylomyelopathy (Wobbler Syndrome)

Authors: Dr. Günter Schwarz & Prof. Dr. Franck Forterre

NOTE: Preoperative radiographs are used to make preliminary selection of spacer and screw sizes. It is advisable but not indispensable to use fluoroscopy while performing the distraction-fusion technique with the C-Lox Cage. Fluoroscopy will provide accurate intraoperative assessment of correct implant and screw size, and of the depth and location of spacer and screw placement.




(1) Place the dog in dorsal recumbency with a fulcrum underneath to support the neck. Avoid overextension. A standard approach to the ventral cervical spine is performed.



(2) The affected disc space is identified and exposed, and a window is cut out of the ventral annulus fibrosus using a beaver blade. The width of this window should be minimally larger than the width of the selected spacer. The vertebral endplates delimit the cranial and caudal borders of the opening. Care must be taken not to penetrate the dorsal part of the annulus during discectomy.



(3) In order to attach the C-LOX Distractor, vertically oriented 2.5 mm holes are drilled into each of the adjacent vertebrae. These holes must be exactly on the midline and should be located in the caudal half of the cranial and in the cranial half of the caudal vertebrae. Penetration depth should at least be 2/3 of the vertebral body depth in order to avoid tilting when distracting the vertebrae. The crista ventralis marks the median plane of the vertebrae but makes it easy to slide off the midline. Preoperative measurements and fluoroscopy will aid in taking care not to enter the spinal canal. Insert the C-LOX Distractor into the predrilled holes and open its jaws to widen the intervertebral disc space.

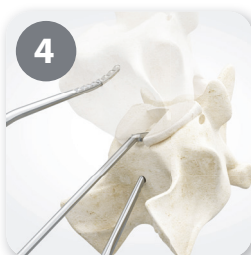

(4) The C-LOX Spinal Disc Broaching Curette is used to carefully remove all remnants of the nucleus pulposus. The dorsal part of the annulus fibrosus can be felt as a more dense structure and should be preserved. The exposed

end plates are freed from as much connective tissue as possible, but subchondral bone must be fully preserved.

If a considerable amount of nucleus pulposus material is located within the spinal canal, this can be attempted to remove with

the help of fine curettes or delicate arthroscopic biopsy forceps.

(5) The suitable size for the C-Lox implant can be estimated pre-operatively on diagnostic imaging. Attach the C-LOX Implantation Placement Rod to the suitable C-LOX Template and insert it in the intervertebral disc space. Ensure the correct fit. Here again, vertical beam fluoroscopy can be very helpful. If insertion requires energetic forcing or seems to achieve only slight distraction, a slimmer or thicker template should be used.



(6) After establishing optimal fit, the Template is replaced by the corresponding C-LOX Cage. Spikes on the cranial and caudal side will ensure a secure seating of the implant, but can make insertion slightly more difficult compared to the C-Lox Template of the same size. Press the cage firmly down into place.



(7) After proper seating of the C-Lox Cage, the C-Lox Distractor can be removed and the Implantation Placement Rod is removed from the cage. The cage will be secured with locking screws.



(8) Select the drill guide according to the screw length intended to use. Attach the Drill Guide to the Implant and use the dedicated 2.0mm drill bit to drill the hole for the Locking Cortical Screw. This drill bit forms a unit with the Drill Guide and will create the hole in the vertebral body to guide the self tapping screw. Screw Lengths are selected pre-operatively from measurements on the preoperative radiographs in precise latero-lateral projection. The four Locking Cortical Screw are inserted using the LeiStar C-Lox Screwdriver. The screws must be well anchored within the vertebral body and must not penetrate the vertebral end plates or the spinal canal. Again, fluoroscopy can be very helpful in



achieving optimal screw placement.

(9) Cancellous bone or bone substitutes can be apposed to the ventral surface of the treated disc space.

Close the soft tissues in a routine manner.

C-LOX Set

Contains:

Sterilization Tray with Lid

1 of each Spinal Cage

1 of each Seizing Cage

2 Placement Rods

Spinal Curette

Screwdriver Shaft

Screwdriver Handle Silicone

1 of each Locking Drill Guide

2 Drills (2.0/2.5mm)


5 of each Screw (10, 12, 14, 16, 18, 20mm)

(30 pcs. total)

134-0100-01

Tray without contents

134-0500-01

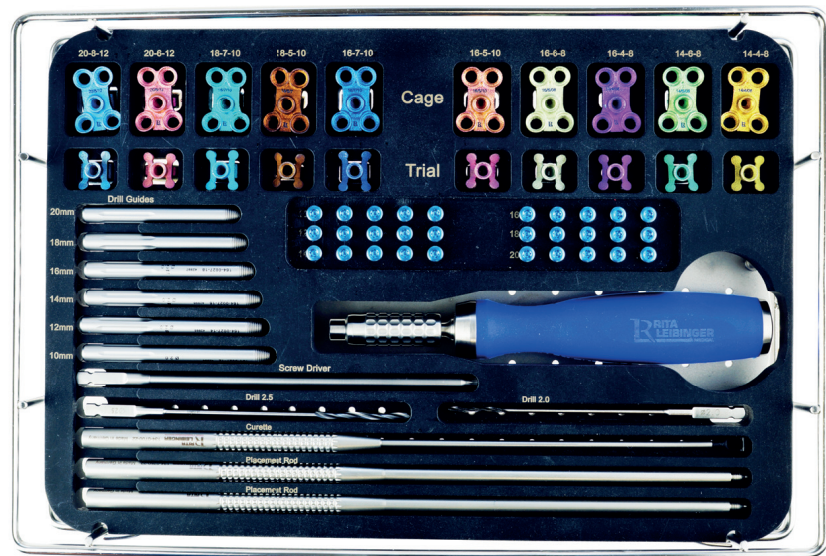

C-LOX Cervical Intervertebral Fusion Implant

C-LOX was developed

by RITA LEIBINGER in cooperation with

Professor Dr. Franck Forterre,

Bern University (Switzerland).



C-LOX Cages

Titanium



14 x 4 x 8 mm



14 x 6 x 8 mm



16 x 4 x 8 mm



16 x 6 x 8 mm



16 x 5 x 10 mm



16 x 7 x 10 mm



18 x 5 x 10 mm



18 x 7 x 10 mm




20 x 6 x 12 mm



20 x 8 x 12 mm

Product Code	Dimensions		
	Depth	Thickness	Width
134-0144-08	14 mm	4 mm	8 mm
134-0146-08	14 mm	6 mm	8 mm
134-0164-08	16 mm	4 mm	8 mm
134-0165-10	16 mm	5 mm	10 mm
134-0166-08	16 mm	6 mm	8 mm
134-0167-10	16 mm	7 mm	10 mm
134-0185-10	18 mm	5 mm	10 mm
134-0187-10	18 mm	7 mm	10 mm
134-0206-12	20 mm	6 mm	12 mm
134-0208-12	20 mm	8 mm	12 mm

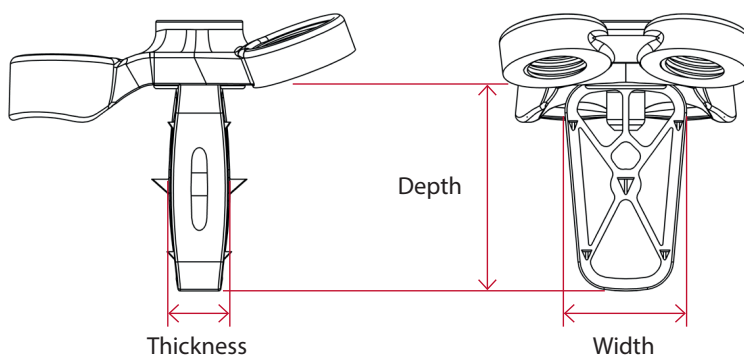

Screws (C-LOX)

C-LOX Templates

Locking Selfdrilling Cortical Screw

Titanium, LeiStar Head, self-tapping with three flute cutting edge




Product Code	Total Length	Thread Length
245-427-08	10 mm	8 mm
245-427-10	12 mm	10 mm
245-427-12	14 mm	12 mm
245-427-14	16 mm	14 mm
245-427-16	18 mm	16 mm
245-427-18	20 mm	18 mm
245-427-20	22 mm	20 mm

C-LOX Templates

Titanium

Product Code	Dimensions		
	Depth	Thickness	Width
134-9144-08	14 mm	4 mm	8 mm
134-9146-08	14 mm	6 mm	8 mm
134-9164-08	16 mm	4 mm	8 mm
134-9165-10	16 mm	5 mm	10 mm
134-9166-08	16 mm	6 mm	8 mm
134-9167-10	16 mm	7 mm	10 mm
134-9185-10	18 mm	5 mm	10 mm
134-9187-10	18 mm	7 mm	10 mm
134-9206-12	20 mm	6 mm	12 mm
134-9208-12	20 mm	8 mm	12 mm

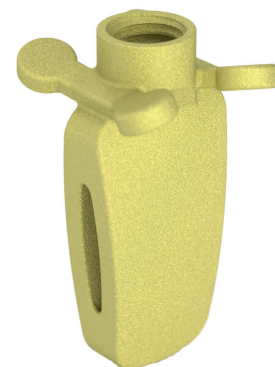

C-LOX Spinal Curette

22cm

134-0700-22



C-LOX Implant Placement Rod

22 cm

134-0750-22



C-LOX Distractor

Product Code	Description
134-0750-16	50mm deep, 160mm length, 90° curved
134-0750-20	50mm deep, 200mm length, 90° curved



Screwdriver Handle

Silicone, AO-Connection
sterilizable up to 134°C / 273°F

247-0102-00



Torque Limiting Screwdriver Handle

Torque 3 Nm
Silicone, AO, sterilizable up to 134°C / 273°F

247-0104-00



C-LOX Screwdriver Shaft 2.7

LeiStar, 135mm length, AO

134-0800-27



C-LOX Locking Drill Guides

Stainless Steel
for 2.0mm drills, for 2.7mm screws



Product Code	Description
164-0027-10	10 mm working end
164-0027-12	12 mm working end
164-0027-14	14 mm working end


Product Code	Description
164-0027-16	16 mm working end
164-0027-18	18 mm working end
164-0027-20	20 mm working end

C-LOX Drill


Product Code	Ø (mm)	Connection
148-0081-20	2.0	AO QC
148-0081-25	2.5	AO QC



LeiCAGE

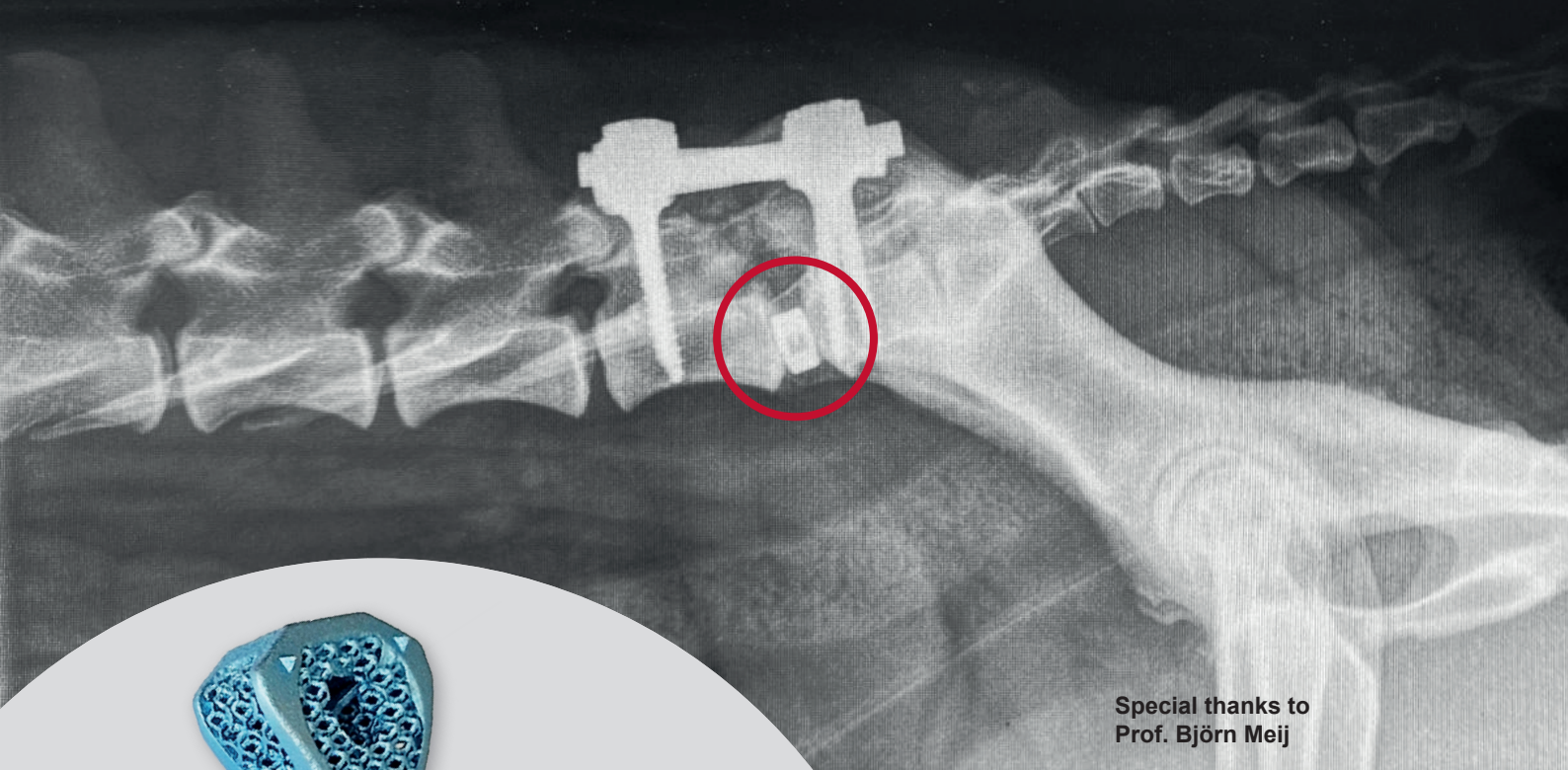

Special thanks to
Prof. Björn Meij

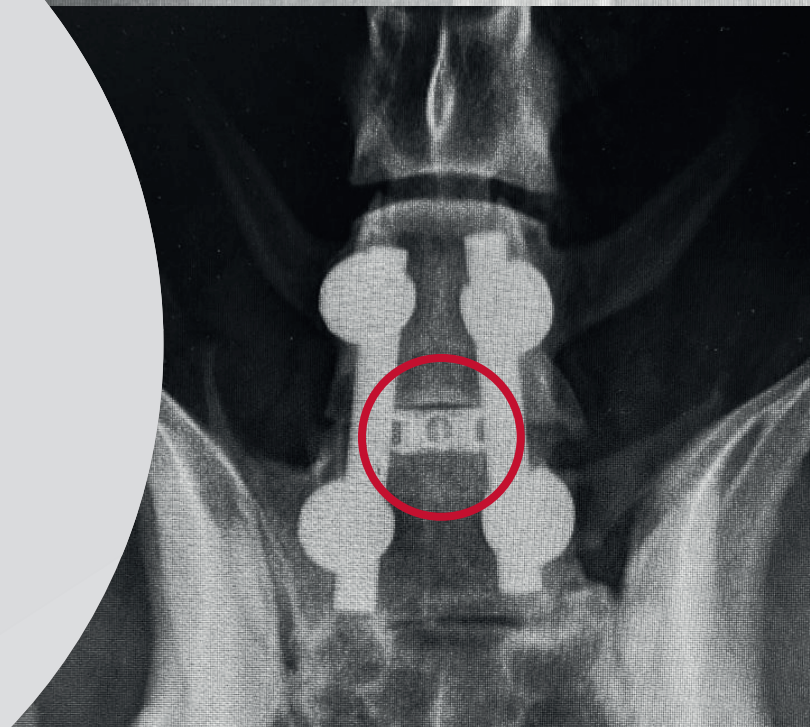

INTERVERTEBRAL CAGES

R RITA
LEIBINGER
MEDICAL

NEW

LeiCAGE Intervertebral Titanium Cages

Restoring Height. Supporting Fusion. Ensuring Stability.

The LeiCAGE Intervertebral Titanium Cage System is designed to support distraction, provide stable intervertebral support, and promote reliable spinal fusion in canine thoracolumbar and lumbosacral procedures. Developed for use in combination with pedicle screw-rod systems such as the new LeiPED, LeiCAGE offers a robust, biologically integrated solution for the stabilization of the vertebral column following decompression or disc removal.

Each LeiCAGE implant is 3D-printed from medical-grade titanium, featuring a porous lattice structure (dodecahedron construction) that encourages rapid bone ingrowth and long-term osseointegration. This design replicates cancellous bone architecture, facilitating direct biological fusion between vertebral bodies and promoting stable load transfer during the healing process.


LeiCAGE implants are available in multiple sizes to accommodate anatomical variations and surgical requirements. Their anatomically contoured geometry ensures optimal endplate contact and uniform force distribution, reducing the risk of subsidence and maximizing stability.

LeiCAGE Intervertebral Cages

Titanium



Product Code	Size	Dimensions (mm)		
		W	D	H
134-0090-07	7 mm	9	7	4.5
134-0105-08	8 mm	10.5	8	5
134-0120-09	9 mm	12	9	5.5
134-0135-10	10 mm	13.5	10	6
134-0150-11	11 mm	15	11	6.5

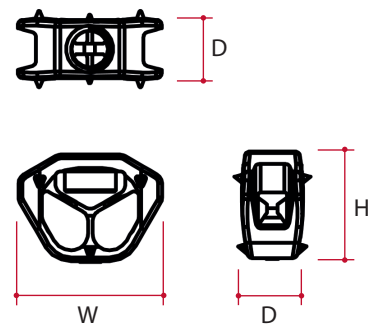

Implant Placement Rod


Stainless Steel




Product Code	Size	Length	Compatibility
134-0751-22	M3	22 cm	For LeiCAGE sizes 7mm, 8mm, 9mm
134-0750-22	M4	22 cm	For LeiCAGE sizes 10mm, 11mm

Lei^{RED}

S/M/L Pedicle Screw System

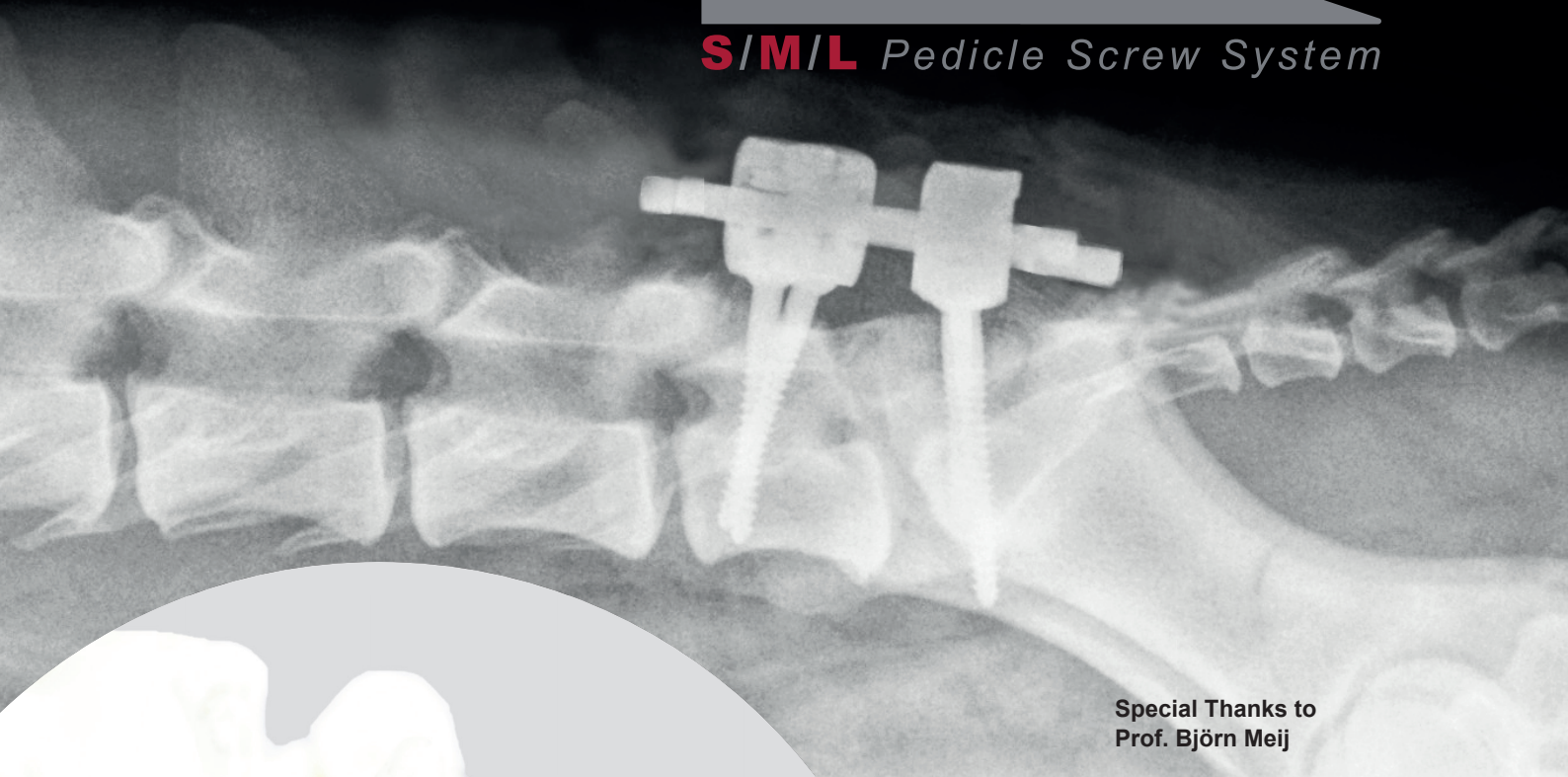

Special Thanks to
Prof. Björn Meij

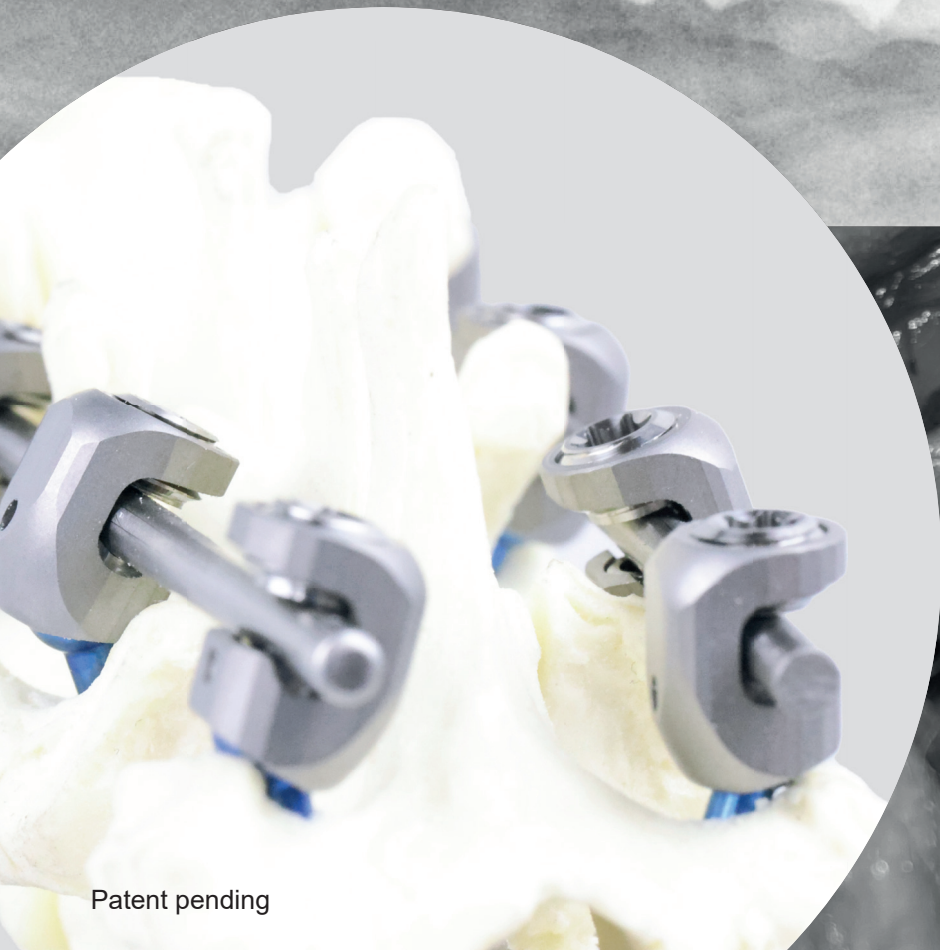

Patent pending



PEDICLE SCREW SYSTEM

RITA
LEIBINGER
MEDICAL

NEW

LeiPED Pedicle Screw System


Advanced Stabilization for the Canine Spine

The **LeiPED Pedicle Screw System** was engineered from the ground up to address the unique anatomical and clinical demands of veterinary spinal surgery. Manufactured from **medical-grade titanium**, it ensures excellent biocompatibility, high mechanical durability, and clear radiographic visibility with minimal imaging artifacts. LeiPED is designed to give veterinary surgeons confidence in treating challenging spinal pathologies such as **degenerative lumbosacral stenosis (DLSS)**, **vertebral instability**, and **spinal trauma**.

GET TO KNOW THE SYSTEM


The LeiPED Pedicle Screw System is a comprehensive fixation solution developed for the unique demands of veterinary spinal surgery. With its innovative design and extensive range of implant options, LeiPED provides surgeons with the versatility, precision, and strength required to treat a wide spectrum of thoracolumbar and lumbosacral spinal conditions.

The system includes tulip screw heads in three size groups – Small (2.0 and 2.4 mm), Medium (2.7 and 3.0 mm), and Large (3.5 and 4.5 mm) – each available in multiple screw lengths to accommodate diverse patient anatomies. Matching rod diameters of 3.0, 4.0, and 5.0 mm ensure stable construct assembly across different case requirements. All screws are supplied with corresponding caps for secure rod fixation. To expand construct versatility, the system also allows the use of cross-connectors, providing additional stability where clinically required.

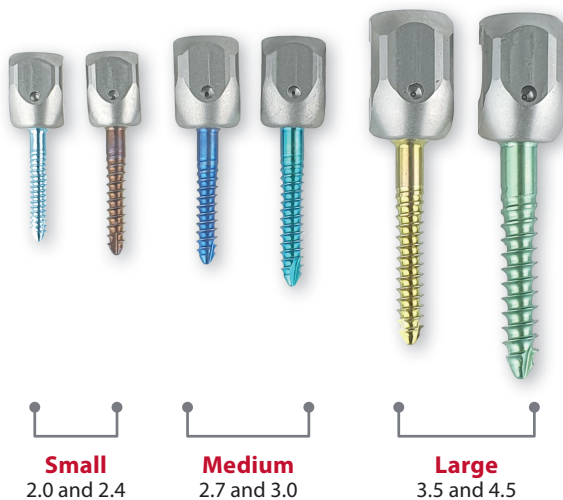


In addition to implants, LeiPED is supported by a complete implant set and dedicated instrument kit, offering surgeons a streamlined and efficient workflow in the operating room. The instruments are designed for intuitive handling, precise screw placement, and reliable rod insertion, ensuring consistent results and surgical confidence.

With its wide selection of implants, compatibility across different spinal pathologies, and ergonomically designed instruments, the LeiPED Pedicle Screw System delivers a full-service solution for advanced veterinary spinal stabilization.

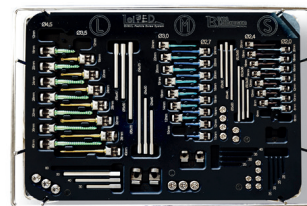

LeiPED Pedicle Screw Sets
Available in **Small**, **Medium**, and **Large**, each set includes pedicle screws in various lengths, matching rods, cross-connectors, and screw caps.

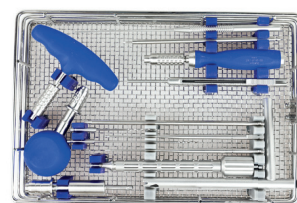

LeiPED Instrument Sets
Also offered in **Small**, **Medium**, and **Large** to match the corresponding implant sets, providing surgeons a size-specific toolkit.

Advantages of LeiPED over Conventional Systems



Side-Loading Rod Placement for Stability

Positioning the rod closer to the bone shortens the lever arm and reduces stress on screw threads and surrounding bone. This biomechanical advantage translates into **lower risk of screw loosening** and **improved force transmission**, ensuring stable vertebral fixation under load.



Asymmetric Screw Head for Anatomical Integration

The asymmetric tulip head adapts to the natural contour of the vertebrae, creating a secure bone-conforming fit. Unlike symmetrical designs, this geometry reduces micro-motion at the bone-implant interface, which in turn minimizes soft tissue irritation.



Variable Polyaxial Angulation for Surgical Flexibility

The screw heads allow controlled polyaxial angulation, giving surgeons the ability to fine-tune screw alignment and rod placement even in anatomically challenging or distorted spinal regions, making implantation safer and more precise.



Compact, Tissue-Sparing Profile

The compact design of the screw head and rod construct facilitates implantation in tight anatomical corridors, reduces implant prominence and protects surrounding musculature and vascular structures from irritation or damage.

High Load-Bearing Capacity with Fatneck Design

The specially engineered Fatneck screw design improves stability under repetitive loading and enhances the construct's resistance to mechanical failure, making LeiPED particularly suited to larger or more active canine patients where higher forces act on the stabilization system.



LeiPED Implants Set - Small 2.0/2.4

LeiStar, Titanium
134-2024-00

Product Code	Description	Qty
134-5220-06	LeiPED Screw 2.0 6mm	2
134-5220-08	LeiPED Screw 2.0 8mm	4
134-5220-10	LeiPED Screw 2.0 10mm	4
134-5220-12	LeiPED Screw 2.0 12mm	4
134-5220-14	LeiPED Screw 2.0 14mm	4
134-5220-16	LeiPED Screw 2.0 16mm	2
134-5220-18	LeiPED Screw 2.0 18mm	2
134-5220-20	LeiPED Screw 2.0 20mm	2
134-5220-22	LeiPED Screw 2.0 22mm	2
134-5224-06	LeiPED Screw 2.4 06mm	2
134-5224-08	LeiPED Screw 2.4 08mm	2
134-5224-10	LeiPED Screw 2.4 10mm	4
134-5224-12	LeiPED Screw 2.4 12mm	4
134-5224-14	LeiPED Screw 2.4 14mm	4
134-5224-16	LeiPED Screw 2.4 16mm	4
134-5224-18	LeiPED Screw 2.4 18mm	2
134-5224-20	LeiPED Screw 2.4 20mm	2
134-5224-22	LeiPED Screw 2.4 22mm	2
134-4225-25	Rod Ø2.5 x 25	4
134-4225-35	Rod Ø2.5 x 35	4
134-4225-50	Rod Ø2.5 x 50	2
134-4225-230	Rod Ø2.5 x 230	1
134-4230-30	Rod Ø3.0 x 30	4
134-4230-40	Rod Ø3.0 x 40	4
134-4230-50	Rod Ø3.0 x 50	2
134-4230-230	Rod Ø3.0 x 230	1
134-4325-20	Cross Connecting Rod Ø2.5 x 20	1
134-4325-25	Cross Connecting Rod Ø2.5 x 25	1
134-4425-00	Cross Connector Tulip - small	4
134-2024-99	Set Screw - small	10
134-2024-01	Tray for Implants Set, small	1

LeiPED Instruments Set - Small


Compatible with the LeiPED Implants Set Small
134-2024-10

Product Code	Description	Qty
134-5520-00	Pedicle finder 2.0	1
134-5524-00	Pedicle finder 2.4	1
134-5507-10	Pedicle screw driver - small	1
134-5508-16	Torque Limiter - small	1
134-5501-10	Couter Torque - small	1
134-5502-10	Rod Holder - medium/small	1
134-5503-10	Set Screw Driver - small	1
134-2024-11	Tray for Instruments Set, small	1

LeiPED Implants Set - Medium 2.7/3.0

LeiStar, Titanium

134-2730-00

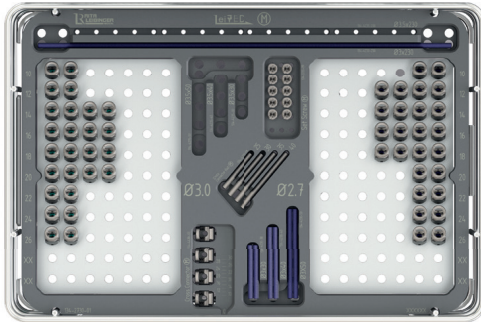


Product Code	Description	Qty
134-5227-10	LeiPED Screw 2.7 10mm	2
134-5227-12	LeiPED Screw 2.7 12mm	4
134-5227-14	LeiPED Screw 2.7 14mm	4
134-5227-16	LeiPED Screw 2.7 16mm	4
134-5227-18	LeiPED Screw 2.7 18mm	4
134-5227-20	LeiPED Screw 2.7 20mm	2
134-5227-22	LeiPED Screw 2.7 22mm	2
134-5227-24	LeiPED Screw 2.7 24mm	2
134-5227-26	LeiPED Screw 2.7 26mm	2
134-5230-10	LeiPED Screw 3.0 10mm	2
134-5230-12	LeiPED Screw 3.0 12mm	2
134-5230-14	LeiPED Screw 3.0 14mm	4
134-5230-16	LeiPED Screw 3.0 16mm	4
134-5230-18	LeiPED Screw 3.0 18mm	4
134-5230-20	LeiPED Screw 3.0 20mm	4
134-5230-22	LeiPED Screw 3.0 22mm	2
134-5230-24	LeiPED Screw 3.0 24mm	2
134-5230-26	LeiPED Screw 3.0 26mm	2
134-4230-30	Rod Ø3.0 x 30	4
134-4230-40	Rod Ø3.0 x 40	4
134-4230-50	Rod Ø3.0 x 50	2
134-4230-230	Rod Ø3.0 x 230	1
134-4240-40	Rod Ø4.0 x 40	4
134-4240-50	Rod Ø4.0 x 50	4
134-4240-70	Rod Ø4.0 x 70	2
134-4240-230	Rod Ø4.0 x 230	1
134-4330-25	Cross Connecting Rod Ø3.0 x 25	1
134-4330-30	Cross Connecting Rod Ø3.0 x 30	1
134-4430-00	Cross Connector Tulip - medium	4
134-2730-99	Set Screw - medium	10
134-2730-01	Tray for Implants Set, medium	1

LeiPED Instruments Set - Medium

Compatible with the LeiPED Implants Set Medium


134-2730-10

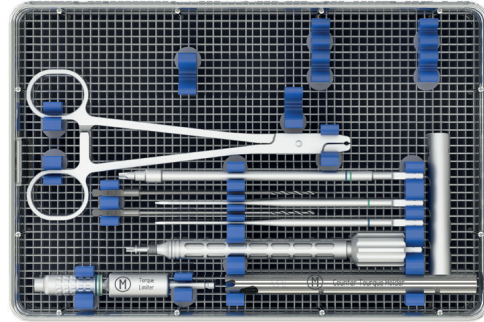

Product Code	Description	Qty
134-5527-00	Pedicle finder 2.7	1
134-5530-00	Pedicle finder 3.0	1
134-5507-20	Pedicle Screw Driver - medium	1
134-5508-30	Torque Limiter - medium	1
134-5501-20	Counter Torque - medium	1
134-5502-10	Rod Holder - medium/ small	1
134-5503-20	Set Screw Driver - medium	1
134-2730-11	Tray for Instruments Set, medium	1

LeiPED Implants Set - Large 3.5/4.5


LeiStar, Titanium
134-3545-00

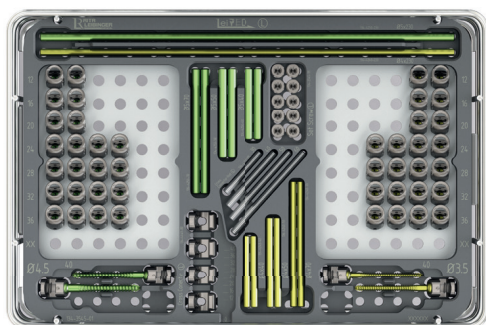

Product Code	Description	Qty
134-5235-12	LeiPED Screw 3.5 12mm	2
134-5235-16	LeiPED Screw 3.5 16mm	4
134-5235-20	LeiPED Screw 3.5 20mm	4
134-5235-24	LeiPED Screw 3.5 24mm	4
134-5235-28	LeiPED Screw 3.5 28mm	4
134-5235-32	LeiPED Screw 3.5 32mm	2
134-5235-36	LeiPED Screw 3.5 36mm	2
134-5235-40	LeiPED Screw 3.5 40mm	2
134-5245-12	LeiPED Screw 4.5 12mm	2
134-5245-16	LeiPED Screw 4.5 16mm	2
134-5245-20	LeiPED Screw 4.5 20mm	4
134-5245-24	LeiPED Screw 4.5 24mm	4
134-5245-28	LeiPED Screw 4.5 28mm	4
134-5245-32	LeiPED Screw 4.5 32mm	4
134-5245-36	LeiPED Screw 4.5 36mm	2
134-5245-40	LeiPED Screw 4.5 40mm	2
134-4240-40	Rod Ø4.0 x 40	4
134-4240-50	Rod Ø4.0 x 50	4
134-4240-70	Rod Ø4.0 x 70	2
134-4240-230	Rod Ø4.0 x 230	1
134-4250-40	Rod Ø5.0 x 40	4
134-4250-50	Rod Ø5.0 x 50	4
134-4250-70	Rod Ø5.0 x 70	2
134-4250-230	Rod Ø5.0 x 230	1
134-4340-30	Cross Connecting Rod Ø4.0 x 30	1
134-4340-40	Cross Connecting Rod Ø4.0 x 40	1
134-4440-00	Cross Connector Tulip - large	4
134-3545-99	Set Screw - large	10
134-3545-01	Tray for Implants Set, large	1

LeiPED Instruments Set - Large


Compatible with the LeiPED Implants Set Large
134-3545-10

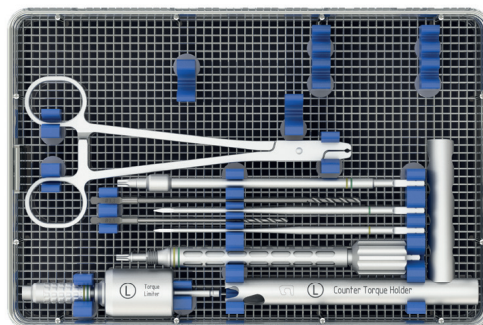

Product Code	Description	Qty
134-5535-00	Pedicle finder 3.5	1
134-5545-00	Pedicle finder 4.5	1
134-5507-30	Pedicle Screw Driver - large	1
134-5508-50	Torque Limiter - large	1
134-5501-30	Counter Torque - large	1
134-5502-30	Rod Holder - large	1
134-5503-30	Set Screw Driver - large	1
134-3545-11	Tray for Instruments Set, large	1

LeiPED Starter Instrument Set

Basic Instrument Set für LeiPED Systems
134-2045-10

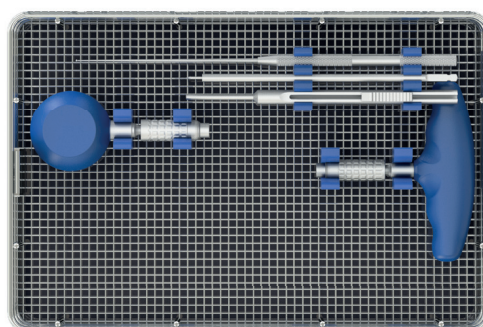

Product Code	Description	Qty
134-2045-11	Tray for Starter Instrument Set	1
164-1520-20	Depth Gauge	1
247-0106-00	Ball Handle AO	1
247-0101-00	T-Handle AO	1
134-5500-20	Pedicle Feeler	1
134-5500-10	Awl with stop	1

ORTHO
4VETS



NEURO
4VETS



Rent your surgery sets



fixed rental fee per surgery + cost of used implants



NEURO
4VETS



POWERTOOLS

hekaDrill

Osteodrive



hekaDrill




ULTIMATE PRECISION

The hekaDrill has been developed with innovation and precision as a guiding principle.

Key capabilities include:

- High speed motor control - up to 85,000 rpm. Allows high torque/low speed control. (min. 200 rpm)
- Advanced motor control ensures comfort for the surgeon during extended use with minimal heat gain.
- Adjustable bur exposure lengths.
- Wireless foot pedal operates using a secure encrypted network, ensuring greater user protection.
- Adjustable acceleration and braking functionality.

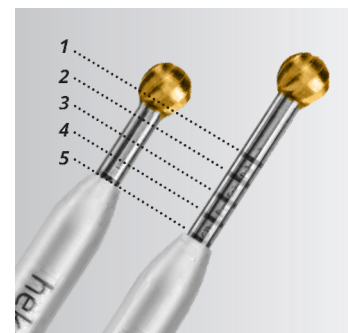

A single rotary bur allows 5 possible exposure positions offering total flexibility and reduced inventory costs.

deSouther
MEDICAL



osteodrive® //

SIMPLICITY, CONTROL & PRECISION



